L I F E I N B A L A N C E

SOEHNLE

VARIO HEAT PILLOW

ACTIVE FIT

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OPERATING INSTRUCTIONS

MODE D'EMPLOI

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BRUGSANVISNING

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IMPORTANT NOTES!

PLEASE STORE THESE INSTRUCTIONS FOR FURTHER USE!

1. Proper use	13
2. Safety information	14
3. Description of the device and operating elements	16
4. Start-up operation	16
5. Temperature control	17
6. Auto Off	17
7. Cleaning and care	18
8. Storage and periodic inspection	19
9. Problems – Causes - Troubleshooting	19
10. Declaration of conformity	20
11. Disposal	20
12. Consumer service	20
13. Guarantee	21
14. Technical information	22

1. Proper use

Thank you for purchasing this device.

With this electronic heat unit you can add effective, targeted heat to your shoulder and neck area. Its use supports blood circulation and therefore the relaxation of tense muscle areas and increases the general well-being.

The device is intended exclusively for domestic use, not for commercial purposes, the use in hospitals or doctor's offices.

In order to be able to utilize all the benefits of your electric warming device, please thoroughly read these instructions and observe the notes on care and safety. You will then be able to enjoy this device for a long time.

Please store these instructions for further use or in case you hand this device over to another person.

Further information can be found at our website at www.soehnle.com.

2. Safety information

Danger

- This device uses supply voltage. Therefore, the same safety regulations apply as for any other electric device.
- Operate the device only with the stated supply voltage.
- Inspect that all components of the device are in a faultless condition prior to each start-up of the device. In the event of a defect or a malfunction, switch the device off immediately and disconnect the power plug from the power network. Repairs must be performed exclusively by authorized customer service technicians. Never open the device. The individual components of this device may not be disassembled, repaired or modified.
- Use the device only in dry rooms and with a dry body – never in the bathtub, shower, above filled sinks, in pools, saunas, or outdoors.
- Never operate the device too close to other heat sources.
- Do not pull the power cord across sharp edges.
- Do not pull on the cord or the device to disconnect the power plug or transport the device.
- Disconnect the power plug if the device is not in use.
- If you dispose of the device at any one point in time, please cut off the cable to make it useless.
- This device may not be used by persons (to include children) with physical, sensory, or mental disabilities or by persons with insufficient experience or knowledge, unless such persons use the device under the supervision of a responsible person or after receiving necessary instructions.
- Please ensure that children do not use the device as a toy.

- Disconnect the device from the power network prior to cleaning it.

During cleaning, ensure that no liquids come into contact with the switch and buttons and that the device is used only in an absolutely dry condition.

- Do not insert any needles or pointy objects into the heat unit.

Warnings

- Please consult a physician prior to using the device if you suffer from an illness (e.g. varicose veins, thrombosis, etc.).
- Never use the device with open wounds, burns, hematoma, swellings, or similar conditions.
- Stop treatment if it feels uncomfortable to you. Heat-sensitive persons in particular should use the device with care and stop the treatment in a timely manner.
- Never use the heat device on animals.
- The use of this device at maximum heat setting for a long period of time can lead to burns of the skin.
- If the device is to be used continuously for a long period of time, the control must be set to the lowest position.

Attention

- In order to ensure the proper functionality of the operating element (thermostat), do not place it on top of or underneath the heat unit.

Do not cover the thermostat with other objects.

- Never use the heat device if it is folded, wrinkled, or pushed together.


3. Device description and operating elements

- Heating element with removable fabric cover
- On/off switch with thermostat and operating light

0 - OFF

1 - 4 Temperature levels

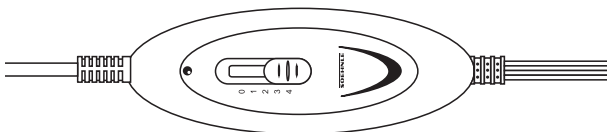4. Start-up operation

Use the Active Fit heat pillow to enjoy comprehensive heat with maximum range of motion.

Its shape is designed for the use in the shoulder and neck area.

Simply place on shoulders and close the hook-and-loop-fasteners.

Important

This heat pillow is equipped with an integrated thermostat, which regulates the temperature automatically and keeps the heat at a constant level. The temperature of the heat pillow corresponds to the maximum permissible value of the applicable IEC regulation regarding user safety. The device generates a mild heat that does not have any harmful effect.

Place the heat pillow in the cover. Ensure that the power cord and the thermostat are outside of the cover. Close the zipper. Before switching the heat pillow on, please ensure that the operating element is in the “0” position. Then insert the power plug into the power outlet.

5. Temperature setting

Use the thermostat to set the desired temperature. The heat pillow switches on automatically when the thermostat is set to a temperature, and the operating light will light up. The heat pillow shuts off and the operating light goes out if the thermostat is set to “0”.

Attention

Please ensure that the heat pillow sits loosely, i.e. some air circulates between the heat pillow and your skin during use.

Accumulation of heat must be prevented. Never use the heat pillow without the cover.

6. Automatic off

The electric heat pillow shuts off automatically after 90 minutes of continuous use. This saves energy and provides additional safety, because the heat pillow is never used without supervision. The operating light will continue to be lit. In order to switch the heat pillow back on, first set the thermostat to “0” and then back to the desired temperature.

After use please disconnect the power plug from the power outlet.

7. Cleaning and care

Cover

The fabric cover of the pillow can be cleaned in the washing machine. Please observe the notes on the labels and the laundry care symbols of the cover.

Disconnect the power plug from the power outlet and remove the heat pillow prior to washing the cover.

- Close the hook-and-loop fasteners for washing.
- Machine or hand-wash at 40°C
- Wash separately or with similar colors.
- Do not bleach
- Iron on low temperature
- No chemical dry-cleaning
- Dry in dryer at low temperature (delicate)

Heat pillow

- Never wash the interior part of the device.
- The electric heat pillow may never be cleaned in a washer nor with chemical agents or at the dry-cleaners.
- The heat pillow cannot be dried in a dryer.
- Never use the heat pillow while it is damp.
- The thermostat cannot come into contact with water under any circumstance. Never hold the heat pillow under running water and do not submerge it into any liquids! The moisture could damage the thermostat.

8. Storage and periodic inspection

- Do not wrinkle the device. Fold the product carefully for storage.
- Store the device at a dry location if it is not in use.
- Let the device cool off prior to folding or rolling it up.
- Inspect the device and power cord regularly for possible signs of wear or damage.
- Do not use chemical moth protection.
- Never place any heavy objects on top of the device as these could damage the electric wires.
- Please return the device to your point of sale in case of any wear or malfunction, where you will receive assistance.

9. Problems - Causes - Troubleshooting

Problem	Cause	Troubleshooting
The device does not heat up.	Device not connected to the power network	Insert power plug into the power outlet and switch device on.
	Device not switched on	Switch on.
	Device defect	During the guarantee period, return the device to your seller who will assist you in receiving a replacement.

10. Declaration of conformity

SOEHNLE herewith declares that this device is in compliance with the basic requirements and all other relevant regulations of applicable European directives. You will find the EU compliance declaration on www.soehnle.com.

11. Disposal of used electronic devices

After its service life, please dispose of this device according to the Directive on Waste Electrical and Electronic Equipment 2002/96/EC (WEEE).

The symbol on the product or its packaging informs you that this product cannot be treated as normal household waste, but must be returned at a collection point for recycling electronic and electrical devices.

12. Customer service

Our customer service will be happy to assist you if you have any questions or comments:

Germany Ph.: (08 00) 5 34 34 34 (free of charge)

International Ph.: +49 26 04 97 70

Monday through Friday 09:00 a.m. to 12:00 noon CET

If you'd like to contact us in writing, please use the contact form at the SOEHNLE homepage www.soehnle.com or use postal mail.

13. Guarantee

SOEHNLE guarantees rectification of faults based on material or production errors, free of charge, either by repair or replacement, for 5 years as of the date of purchase. Please retain the receipt of purchase and guarantee slip. In case of a guarantee claim, please return the device, guarantee slip and receipt of purchase to your seller.

Guarantee slip

Name, first name _____

Street _____

Postal code, city/town _____

Purchase date _____

Reason for claim _____

Seller _____

14. Technical information

Dimensions	~ 47 x 31 cm
Rated power	45 W
Rated voltage	220 - 240 V ~ 50 Hz

The device is equipped with double protective insulation and corresponds to protection class II.

All SOEHNLE heat products are marked with the Eco-Tex Certificate 100, which is globally recognized and certifies that none of the textile components contain any substances that are harmful to your health.

The SOEHNLE Safety Concept guarantees the consumer the maximum level of safety. Exclusive manufacture, additional test criteria, and high SOEHNLE safety standards, e.g. switch monitoring, heating element insulation, auto-off timer, mains-in safety, and thermal safety with auto-off function, convey a good feeling all around.

With this certification mark, the products as well as the production processes are subject to periodic control by KEMA. KEMA monitors conformity continuously. The tests are in compliance with European Standards. KEMA-KEUR is recognized in all European countries.

Type inspection pursuant to European Safety Standards. Complete inspection of the labeling and warnings, construction, mechanical tests, wiring, protection against electric shock, moisture and dust, durability and thermal tests as well as the durability against heat, fire, and leakage current.

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Leifheit AG